

DOCTOR INSTRUCTIONS

Thank you for purchasing the Vision Assessment Corporation Frog Polarized Variable Vectograph Vision Therapy System, P/N 1058.



PURPOSE

A three-dimensional picture is used to strengthen the binocularity system and provide base-in and/or base-out training.

FAMILIARIZE YOURSELF WITH THE SYSTEM

- Polarized Variable Vectograph Vision Therapy System consists of:

1. 1 Guide



2. 2 Vectographic Panels

- Panel 1 (Blue Bar at Bottom)



- Panel 2 (Numbers/Letters at Bottom)



3. 1 Pair Standard Polarized Viewers



(NOT TO BE USED AS SUNGLASSES)

4. 2 Instruction Manuals

- Doctor Instruction Manual
- Patient Instruction Manual

This exercise works best if practiced daily as suggested. A Patient Instruction Manual has been included so the Vectograph Vision Therapy System can be sent home with the patient to allow for greater success through continued at-home usage and practice.

5. 1 Patient Vision Therapy Record Form with Pen

TESTING CONDITIONS

- Well-lit, glare-free area
- If reflections or glare on the Vectograph can be seen, try tilting it or choose another testing location.

ADMINISTRATION

1. Place the polarized viewers on the patient.
2. **PLEASE NOTE:** DOCTOR SHOULD ADVISE WHETHER OR NOT POLARIZED VIEWERS SHOULD BE WORN OVER PATIENT'S PRESCRIPTION GLASSES.
3. Hold the Vectograph approximately 16 inches (40 cm) in front of the patient.
4. Slide the 2 panels in the guide until the 0 appears in the opening on the blue bar at the bottom of the test.
5. Now ask the patient to look at the box containing the + located in the upper right-hand side of the Vectograph.
6. The patient should be able to see both the L & R (Left & Right) below the +. If the patient is unable to see both letters, instruct him/her to cover the opposite eye for the side he/she cannot see. For example, if the patient cannot see the R (Right), have him/her cover his/her left eye. Have the patient focus with the one eye until he/she is able to see the letter that was not previously seen when viewing the Vectograph binocularly. Then again have the patient view the Vectograph binocularly to ensure that now he/she can see both the L & R. Once both the L & R can be seen, the patient is ready to use the Vectograph.
7. Slowly slide the panels apart and ask the patient to try to maintain a single Frog image. *As the panels are separated and the blue bar moves along the letters from A to P, the Frog image should get larger and float behind the Vectograph (divergence). As the panels are separated and the blue bar moves along the numbers from 1 – 24, the Frog image should get smaller and float in front of the Vectograph (convergence).*
8. Have the patient notify you when he/she is no longer able to fuse the Frog image. (This is the Breaking Point).
9. Note the number/letter in the space on the blue bar at the bottom of the Vectograph. Refer to "Scoring" section of this manual.
10. Next, as you slowly slide the panels back together, ask the patient to notify you when he/she is able to see the Frog image as a single clear image again. (This is the Recovery Point).
11. Note the number/letter in the in the space on the blue bar at the bottom of the Vectograph. Refer to "Scoring" section of this manual.
12. Repeat steps 7 – 11 until patient can achieve the doctor's recommended goal.

SUPPLEMENTAL EXERCISE

13. The Frog Vectograph can also be used in conjunction with the Chain Variable Vectograph or Gem Variable Vectograph to provide additional exercises.

Gem Variable Vectograph



(P/N 1060)

Chain Variable Vectograph



(P/N 1056)

14. Hold the Frog Vectograph behind the Gem/Chain Vectograph and ask the patient to look through the Gem/Chain image while viewing the Frog image behind and centered inside of the Gem or Chain.
15. Repeat steps 7 –11 using the Gem/Chain Vectograph while holding the Frog Vectograph stationary behind the Gem/Chain as the patient continues to view the Frog image through the Gem/Chain.
16. The use of the additional Vectograph increases the depth of the three-dimensional images to further strengthen the binocularity system.

SCORING

- Each letter on the bottom blue bar represents one diopter (Base-In) (Divergence / Relaxing).
- Each number 1-10 on the bottom blue bar represents one diopter. (Base-Out) (Convergence / Crossing).
- Each even number 10-24 represents two diopters. (Base-Out) (Convergence / Crossing).
- 40 diopter range of separation available.
- Vectograph has a range of 5300 sec of arc.
- Frog image has a range of 2500 sec of arc.
- Refer to the Patient Vision Therapy Record Form.
- Record the date.
- Above the diagonal line in the row labeled “Relaxing” record the Breaking Point for the highest letter that the patient is able to reach.
- In the same row below the diagonal line record the Recovery Point for the highest letter at which the patient is able to fuse the images as a single clear image again.
- Repeat the previous 2 steps for the row labeled “Crossing” for recording the Breaking Point and Recovery Point for the highest numbers that the patient can reach on the blue bar at the bottom of the Vectograph.


CARE/HANDLING & STORAGE

- Clean vectographic panels and guide with a soft, damp, lint-free cloth. Dampen cloth using glass cleaner or mild detergent/water.
- **CAUTION: DO NOT IMMERSE THE VECTOGRAPHIC PANELS IN WATER. DO NOT SPRAY CLEANER DIRECTLY ONTO PANELS.**
- ☀ ☂ Store Vectograph in a dry place away from direct sunlight.
- Clean polarized viewers using lens cleaner and soft, lint-free cloth.
- If panels are removed from guide during cleaning, replace the panels in the guide placing the panel with the blue bar on top of the panel with the numbers/letters and ensuring that the plastic portion of the guide is behind the panels.

WARRANTY

- 1 year manufacturer warranty from date of purchase.

Notice to User/Patient: *Any serious incident that has occurred in relation to this device should be reported to the manufacturer and to the competent authority of the Member State in which the user and/or patient is established.*

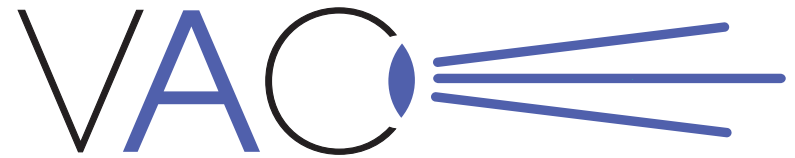

Vision
Assessment
 Corporation TM

The Frog

Polarized Variable Vectograph

P/N 1058
DOCTOR INSTRUCTIONS

Vision Assessment Corporation
5400 Newport Drive, Suite 3
Rolling Meadows, Illinois 60008 USA
Phone: 1 847 239 5889
Email: [sales @VisionAssessment.com](mailto:sales@VisionAssessment.com)
Web: www.VisionAssessment.com



MDSS GmbH,
 Schiffgraben 41, Hannover 30175, Germany

Manufactured in USA by Vision Assessment Corporation © 2013