

DOCTOR INSTRUCTIONS

Thank you for purchasing the Vision Assessment Corporation Chessboard Polarized Variable Vectograph (P/N 1042).



PURPOSE

A three-dimensional picture is used to strengthen the binocularity system and provide base-in and/or base-out training.

FAMILIARIZE YOURSELF WITH THE VECTOGRAPH

- Polarized Variable Vectograph Vision Therapy System consists of:

1. 1 Guide



2. 2 Vectographic Panels

- Panel 1 (Blue Bar at Bottom)



- Panel 2 (Numbers/Letters at Bottom)



3. 1 Pair Standard Polarized Viewers



(NOT TO BE USED AS SUNGLASSES)

4. 2 Instruction Manuals

- Doctor Instruction Manual
- Patient Instruction Manual

This exercise works best if practiced daily as suggested. A Patient Instruction Manual has been included so the Vectograph Vision Therapy System can be sent home with the patient to allow for greater success through continued at-home usage and practice.

5. 1 Patient Vision Therapy Record Form with Pen

TESTING CONDITIONS

- Well-lit, glare-free area
- If reflections or glare on the Vectograph can be seen, try tilting the test or choose another testing location

ADMINISTRATION

1. Place the polarized viewers on the patient.

PLEASE NOTE: DOCTOR SHOULD ADVISE WHETHER OR NOT POLARIZED VIEWERS SHOULD BE WORN OVER PATIENT'S PRESCRIPTION GLASSES.


2. Hold the Vectograph approximately 16 inches (40 cm) in front of the patient.
3. Slide the 2 panels in the guide until the 0 appears in the opening on the blue bar at the bottom of the test.
4. Begin by performing a suppression check. This Vectograph provides two options, and either can be used to ensure the patient is viewing the Vectograph binocularly.
 - a. The patient should be able to see both the L & R (Left & Right) below the +. If the patient is unable to see both letters, instruct him/her to cover the opposite eye for the side he/she cannot see. For example, if the patient cannot see the R (Right), have him/her cover his/her left eye. Have the patient focus with the one eye until he/she is able to see the letter that was not previously seen when viewing the Vectograph binocularly. Then again have the patient view the Vectograph binocularly to ensure that now he/she can see both the L & R. Once both the L & R can be seen, the patient is ready to use the Vectograph.
 - b. Ask the patient to point to/indicate the number of captured Chess Pieces on both the Left of the Chess board and to the Right of the board. There are 6 captured Pieces on the Left and 5 captured Pieces on the Right. If the Patient is unable to see all captured Pieces on each side, have the patient cover the opposite eye for the side which all captured Pieces are not seen. For Example, if the patient only sees 4 Pieces on the Left, the patient is suppressing vision in the Left Eye. Have the patient cover his/her Right Eye until all 5 captured Pieces on the Left can be seen, and then again have the patient view the Vectograph binocularly to ensure all 11 captured Pieces can be seen. When all 11 captured pieces can be seen, the patient is ready to use the Vectograph.
5. The individual Chess Pieces within the Vectograph are arranged at different levels of depth with varying seconds of arc. Ask the patient to point to or name each Chess Piece and indicate whether it appears to be floating away from them (BELOW the Vectograph), towards them (ABOVE the Vectograph) or if it does not appear to float (EVEN with the Vectograph). Have the patient indicate which Chess Pieces appear to be floating away (BELOW) the most and which appear to be floating towards them (ABOVE) the most.
6. Refer to the "Answer Key" section of this manual for the seconds of arc and direction and depth of float for each of the images.
7. Next, while slowly sliding the panels apart ask the patient to try to maintain a single image while repeating Step 5.
8. Have the patient notify you when the images double or when he/she is no longer able to fuse the Chess Piece images. (This is the Breaking Point.)
9. Note the number/letter in the space on the blue bar at the bottom of the Vectograph. Refer to "Scoring" section of this manual.
10. Next, as you slowly slide the panels back together, ask the patient to notify you when he/she is able to see the doubled images as a single clear image again. (This is the Recovery Point.)
11. Note the number/letter in the in the space on the blue bar at the bottom of the Vectograph. Refer to "Scoring" section of this manual.
12. Repeat steps 5 – 11 until patient can achieve the doctor's recommended goal.
13. Follow any variation on these exercises as recommended by the doctor.

SCORING

- Each letter A-P on the bottom blue bar represents one diopter (Base-In) (Divergence / Relaxing).
- Each number 1-10 on the bottom blue bar represents one diopter. (Base-Out) (Convergence / Crossing).
- Each number 10-24 represents two diopters. (Base-Out) (Convergence / Crossing).
- 40 diopter range of separation available.

- Refer to the Patient Vision Therapy Record Form.
- Record the date.
- Above the diagonal line in the row labeled “Relaxing” record the Breaking Point for the highest letter on the blue bar that the patient is able to reach.
- In the same box, below the diagonal line, record the Recovery Point for the highest letter from the blue bar at which the patient is able to fuse the images as a single clear image again.
- Repeat the previous 2 steps for the row labeled “Crossing” for recording the Breaking Point and Recovery Point for the highest numbers that the patient can reach on the blue bar at the bottom of the Vectograph.

ANSWER KEY

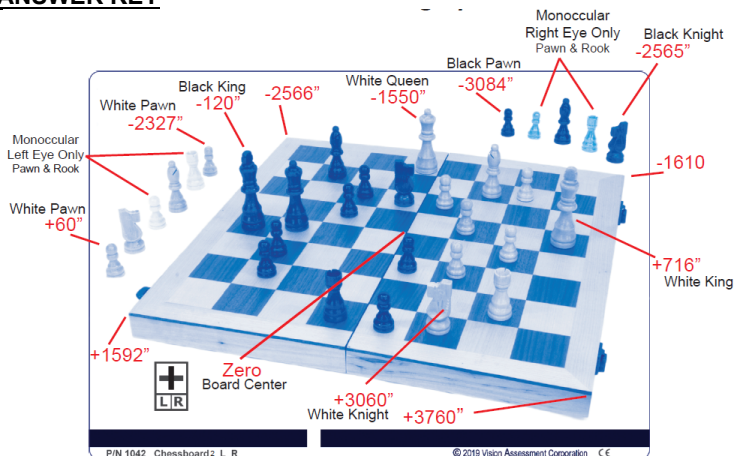

Depth Levels (Seconds of Arc) / Image Orientation to Vectograph: (Clockwise Left to Right)

- | | |
|---|--|
| White Pawn, <i>Front Left</i> : ABOVE | Black Knight, <i>Back Right</i> : BELOW |
| White Pawn, <i>Back Left</i> : BELOW | White King, <i>Right Center</i> : ABOVE |
| Black King, <i>Left Center</i> : BELOW | Board, <i>Front</i> : ABOVE |
| White Queen, <i>Back Center</i> : BELOW | White Knight, <i>Front Right</i> : ABOVE |
| Black Pawn, <i>Back Right</i> : BELOW | Board, <i>Center</i> : EVEN |


CARE/HANDLING & STORAGE

- Clean vectographic panels and guide with a soft, damp, lint-free cloth. Dampen cloth using glass cleaner or mild detergent/water.
- **CAUTION: DO NOT IMMERSE THE VECTOGRAPHIC PANELS IN WATER. DO NOT SPRAY CLEANER DIRECTLY ONTO PANELS.**
- ☀️☔ Store Vectograph in a dry place away from direct sunlight.
- Clean polarized viewers using lens cleaner and soft, lint-free cloth.
- If panels are removed from guide during cleaning, replace the panels in the guide, placing the panel with the blue bar on top of the panel with the numbers/letters. and ensuring that the plastic portion of the guide is behind the panels.

WARRANTY

- 1-year manufacturer warranty from date of purchase

Notice to User/Patient: Any serious incident that has occurred in relation to this device should be reported to the manufacturer and to the competent authority of the Member State in which the user and/or patient is established.

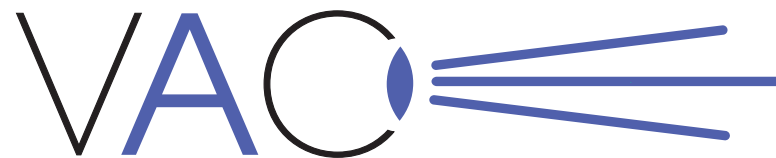

Vision
Assessment
Corporation™

Chessboard

Polarized Variable Vectograph

P/N 1042

DOCTOR INSTRUCTIONS

Vision Assessment Corporation
5400 Newport Drive, Suite 3
Rolling Meadows, Illinois 60008 USA
Phone: 1 847 239 5889
Email: sales@VisionAssessment.com
Web: www.VisionAssessment.com



MDSS GmbH,
Schiffgraben 41, Hannover 30175, Germany

Manufactured in USA by Vision Assessment Corporation © 2019