

## PATIENT INSTRUCTIONS

Thank you for purchasing the Vision Assessment Corporation Pool Table Polarized Variable Vectograph (P/N 1049).



### PURPOSE

A three-dimensional picture is used to strengthen the binocularity system and provide base-in and/or base-out training.

### FAMILIARIZE YOURSELF WITH THE VECTOGRAPH

- Polarized Variable Vectograph Vision Therapy System consists of:

1. 1 Therapy Binder



2. 1 Guide



3. 2 Vectographic Panels

- Panel 1 (Blue Bar at Bottom)



- Panel 2 (Numbers/Letters at Bottom)



4. 1 Pair Standard Polarized Viewers



(NOT TO BE USED AS SUNGLASSES)

5. 2 Instruction Manuals

- Doctor Instruction Manual
- Patient Instruction Manual

*This exercise works best if practiced daily as suggested. A Patient Instruction Manual has been included so the Vectograph Vision Therapy System can be sent home with the patient to allow for greater success through continued at-home usage and practice.*

6. 1 Patient Vision Therapy Record Form with Pen

### TESTING CONDITIONS

- Well-lit, glare-free area
- If reflections or glare on the Vectograph can be seen, try tilting the test or choose another testing location

## ADMINISTRATION

1. Place the polarized viewers on the patient.


**PLEASE NOTE: DOCTOR SHOULD ADVISE WHETHER OR NOT POLARIZED VIEWERS SHOULD BE WORN OVER PATIENT'S PRESCRIPTION GLASSES.**

2. Hold the Vectograph approximately 16 inches (40 cm) in front of the patient.
3. Slide the 2 panels in the guide until the 0 appears in the opening on the blue bar at the bottom of the test.
4. Begin by performing a suppression check. The patient should be able to see both the L & R (Left & Right) below the +. If the patient is unable to see both letters, instruct him/her to cover the opposite eye for the side he/she cannot see. For example, if the patient cannot see the R (Right), have him/her cover his/her left eye. Have the patient focus with the one eye until he/she is able to see the letter that was not previously seen when viewing the Vectograph binocularly. Then again have the patient view the Vectograph binocularly to ensure that now he/she can see both the L & R. Once both the L & R can be seen, the patient is ready to use the Vectograph.
5. Each of the Balls and portions of the Pool Table within the Vectograph are arranged at different levels of depth with varying seconds of arc. Ask the patient to identify each Ball and indicate if it appears to be floating away from them (BELOW the Vectograph), towards them (ABOVE the Vectograph) or if it does not appear to float (EVEN with the Vectograph).
6. Refer to the "Answer Key" section of this manual for the seconds of arc and direction and depth of float for each of the images. Ask the patient to compare various Balls or portions of the Pool Table and indicate which object is closer and which is further away.
7. Next, while slowly sliding the panels apart ask the patient to try to maintain a single image while repeating Steps 5 & 6.
8. Have the patient indicate when the images double or when he/she is no longer able to fuse the Balls and Pool Table image. (This is the Breaking Point.)
9. Note the number/letter in the space on the blue bar at the bottom of the Vectograph. Refer to "Scoring" section of this manual.
10. Next, as you slowly slide the panels back together, ask the patient to indicate when he/she is able to see the doubled images as a single clear image again. (This is the Recovery Point.)
11. Note the number/letter in the in the space on the blue bar at the bottom of the Vectograph. Refer to "Scoring" section of this manual.
12. Repeat steps 5 – 12 until patient can achieve the doctor's recommended goal.
13. Follow any variation on these exercises as recommended by the doctor.

### SCORING

- Each letter A-P on the bottom blue bar represents one diopter (Base-In) (Divergence / Relaxing).
- Each number 1-10 on the bottom blue bar represents one diopter. (Base-Out) (Convergence / Crossing).
- Each number 10-24 represents two diopters. (Base-Out) (Convergence / Crossing).
- 40 diopter range of separation available.
- Refer to the Patient Vision Therapy Record Form and record the date
- Above the diagonal line in the row labeled "Relaxing" record the Breaking Point for the highest letter on the blue bar that the patient is able to reach.
- In the same box, below the diagonal line, record the Recovery Point for the highest letter from the blue bar at which the patient is able to fuse the images as a single clear image again.
- Repeat the previous 2 steps for the row labeled "Crossing" and record the Breaking Point and Recovery Point for the highest number that the patient can reach on the blue bar at the bottom of the Vectograph.

## ANSWER KEY

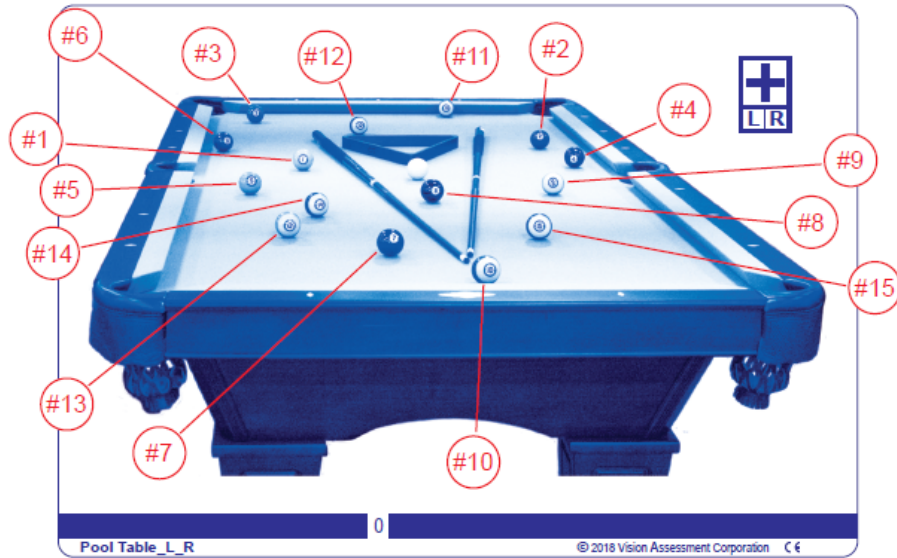

Depth Levels (Seconds of Arc) / Image Orientation to Vectograph:  
(Clockwise Left to Right) Total Stereopsis Depth: 14,365"

Front Rail: 6722" / ABOVE

#10: 4753" / ABOVE

#7: 2684" / ABOVE

#13: 1273" / ABOVE

#14: 139" / BELOW

#5: 1809" / BELOW

#1: 3422" / BELOW

#6: 4934" / BELOW

#3: 6822" / BELOW

Back Rail: 7643" / BELOW

#12: 5907" / BELOW

Cue Ball: 2863 / BELOW

#11: 7134" / BELOW

#2: 1691" / BELOW

#4: 1451" ABOVE

#9: 7134" / BELOW

#8: 1037" / BELOW

#15: 1451" / ABOVE

### CARE/HANDLING & STORAGE

- Clean vectographic panels and guide with a soft, damp, lint-free cloth. Dampen cloth using glass cleaner or mild detergent/water.
- **CAUTION: DO NOT IMMERSE THE VECTOGRAPHIC PANELS IN WATER. DO NOT SPRAY CLEANER DIRECTLY ONTO PANELS.**
- ☀️ ☔ Store Vectograph in a dry place away from direct sunlight.
- Clean polarized viewers using lens cleaner and soft, lint-free cloth.
- If panels are removed from guide during cleaning, replace the panels in the guide, placing the panel with the blue bar on top of the panel with the numbers/letters. and ensuring that the plastic portion of the guide is behind the panels.

### WARRANTY

- 1 year manufacturer warranty from date of purchase



## Pool Table

Polarized Variable Vectograph

**P/N 1049**

**PATIENT INSTRUCTIONS**

**Vision Assessment Corporation**

**2675 Coyle Avenue**

**Elk Grove Village, Illinois 60007 USA**

**Phone: 1 847 239 5889**

**Toll Free USA: 1 866 887 9692**

**Email: [sales@VisionAssessment.com](mailto:sales@VisionAssessment.com)**

**Web: [www.VisionAssessment.com](http://www.VisionAssessment.com)**



MDSS GmbH,  
Schiffgraben 41, Hannover 30175, Germany

Manufactured in USA by Vision Assessment Corporation © 2018